

CoinTip™ snare



 **STERIS**

Endoscopy

The logo for CoinTip SNARE features the words "CoinTip" in a large, bold, blue sans-serif font, with "SNARE" in a smaller, blue sans-serif font below it. A blue line representing a snare wire curves over the text, ending in a loop that resembles a coin tip. The background is white with a green header at the top and several blue curved lines on the right side.

CoinTip

SNARE

STERIS Endoscopy is continuously innovating to deliver new tools and technologies in the polypectomy and tissue acquisition space.

Our newest addition, the CoinTip snare, is a truly differentiated solution for complex polypectomy procedures including endoscopic mucosal resection (EMR).



PRODUCT DETAILS:

Not all snares offer the same recruitment, resection, and performance. The innovative CoinTip snare supports hot and cold resection with very unique design features.

Coined Distal Tip / Snare Loop Wire:

The CoinTip snare leverages the benefits of different wire diameters to support precise hot and cold cutting capabilities.

- Distal tip flattened to 0.010" which is achieved by the manufacturing process of "coining" (*patent pending feature*)
- 0.016" snare loop wire reinforces grip of tissue to capture target
- Snare loop maintains its shape and integrity after multiple resections

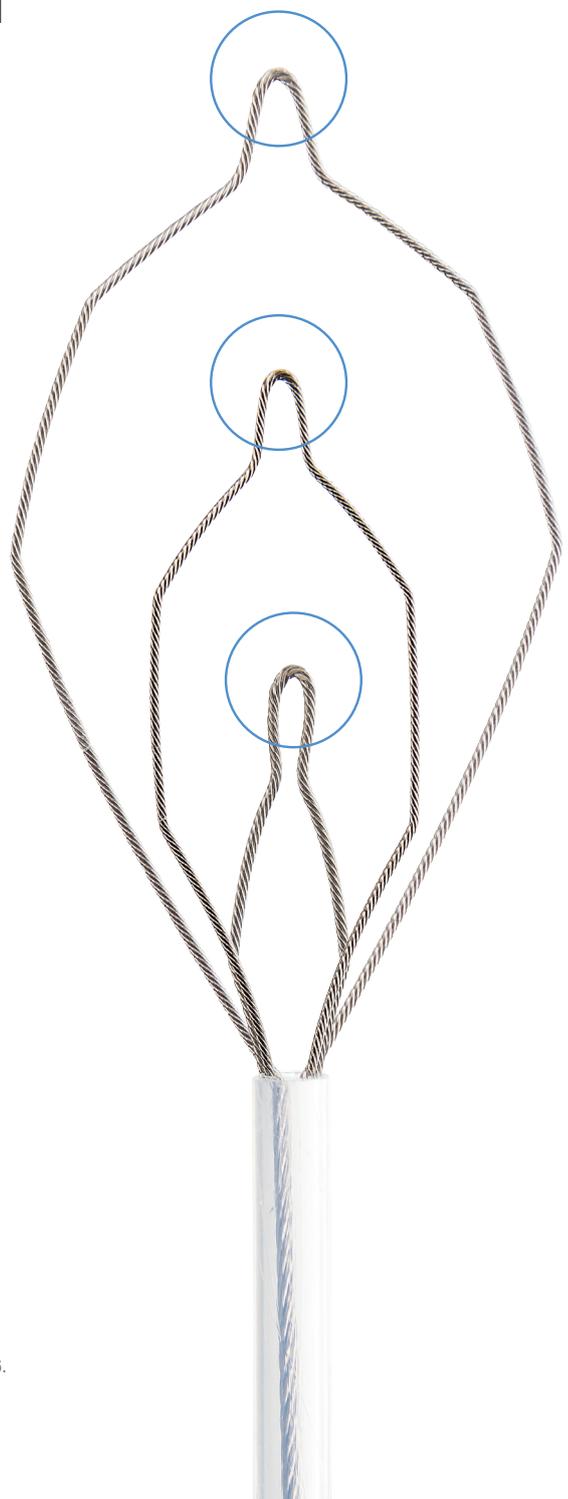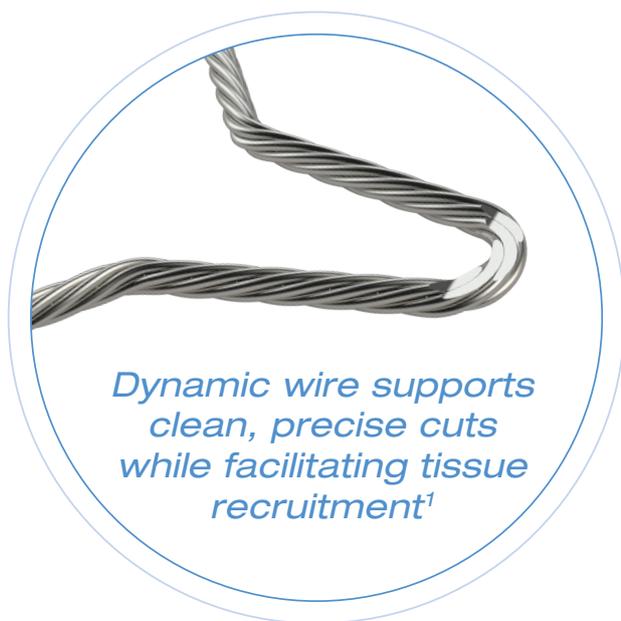


¹Din S, Ball A. "Cold snare polypectomy: Does snare type influence outcomes?" Digestive Endoscopy (2015): 1-6.

Double Catheter:

- Catheter includes an internal strain relief for added support during resection to reduce the chance of buckling during resection
- Traditional snare catheters often buckle resulting in incomplete resection¹


¹Hewett D. "Cold snare polypectomy: optimizing technique and technology." *Gastrointestinal Endoscopy* 82 (2015): 693-696.

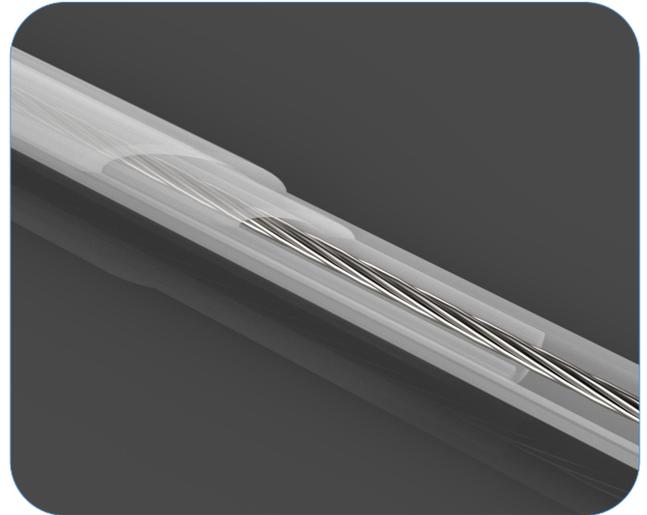

Multiple Stages / Memory Points:


The CoinTip snare has dedicated memory points to support EMR and a variety of resection needs, including multiple resections or various types and sizes of polyps.

- **20mm Deltoid Snare** – ability to encompass majority of lesions with control
- **10mm Hexagonal Snare** – offers stiffness and precision for additional control
- **6mm Diamond Snare** – specific placement for diminutive polyps and cleaning up resection site margins

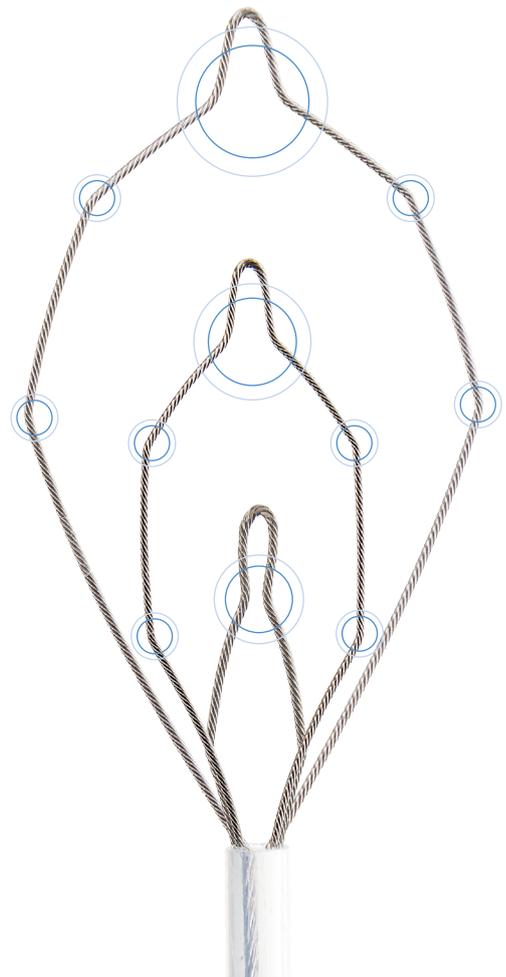


EMR CASES:

The unique features of the CoinTip snare make it the ideal device for physicians performing challenging polypectomies. The CoinTip snare has dual indications for use with and without cautery providing technique versatility needed to resect a variety of polyps.

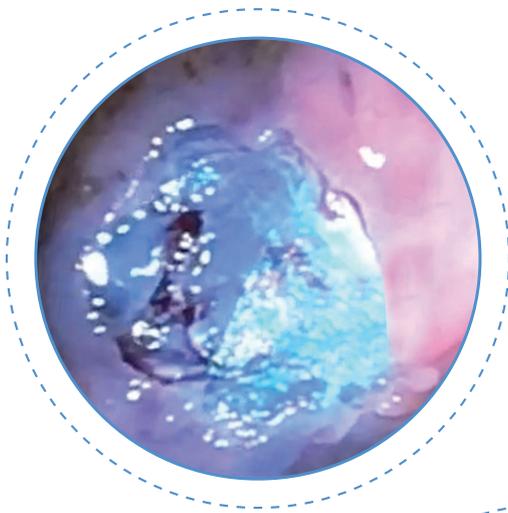

*Resection site after hot EMR
with minimal cautery effect*

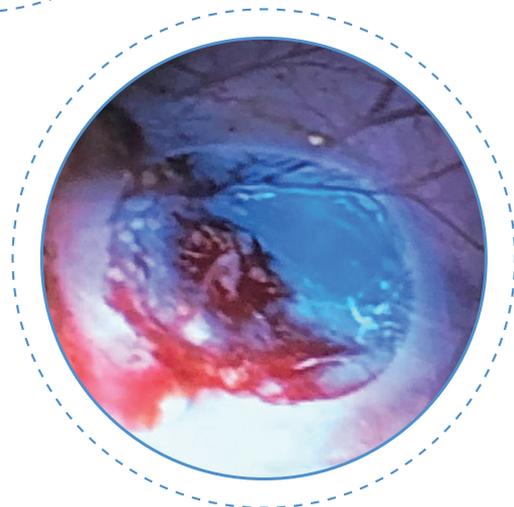

*Resection site after cold
EMR en bloc resection*



ORDERING INFORMATION:

For more information on the CoinTip snare, contact STERIS Endoscopy today or visit www.steris.com.

Product Number	Description	Sheath Diameter (mm)	Sheath Length (cm)	Approx. Snare Size (W x L)	Snare Wire Diameter	Units/Box
BX00711110	CoinTip snare	2.5	230	Deltoid Shape: 20 x 32mm Hexagonal Shape: 10 x 27mm Diamond Shape: 6 x 19mm	0.41mm (0.016") Tip: 0.25mm (0.010")	10

Distributed by:



©2020 STERIS Endoscopy. All Rights Reserved.
Unless otherwise indicated, all marks denoted with ® or ™ are registered with the U.S. Patent and Trademark Office, or are trademarks owned by STERIS or its affiliates.

Follow us:    

STERIS
5960 Heisley Road
Mentor, OH 44060-1834 – USA
Customer Service: 800-548-4873
Fax: 440-639-4495
www.steris.com